



Secret™ PRO

Progressive Dermal
Remodeling with
Defined Results

CUTERA®
FACE + BODY AESTHETIC SOLUTIONS

Secret™ PRO


Take Anti-aging to New Depths

Secret PRO uniquely combines two proven technologies, radio frequency microneedling and fractional CO₂, into a single device for skin resurfacing and deep dermal remodeling.

Treat the skin from the epidermis to the deeper dermal layers, independently or together, with Secret PRO's signature treatment options. Customize treatments based on each individual patient's needs, skin type, and downtime preferences to improve patient experience, results, and satisfaction.



2 in 1

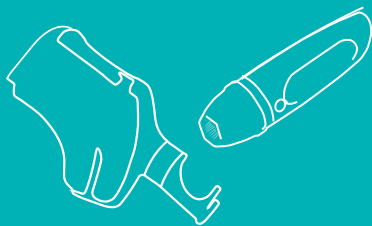

RF microneedling and fractional CO₂ laser for comprehensive skin resurfacing and dermal remodeling.

Customizable



Set needle penetration depth, RF energy, fractional ablation, scan pattern, and scan mode.

Excellent ROI




Treat all skin types with semi-insulated needles and offer CO₂ and non-insulated needles to skin types 1-3. Most treatments are quick and delegable with high patient throughput and satisfaction.

RF + CO₂: A Dynamic Duo

2 MHz Bipolar Fractional Radio Frequency Microneedling

Secret PRO's advanced bipolar RF energy denatures the tissue surrounding each microneedle stimulating collagen, hyaluronic acid production, and dermal remodeling through hemostasis and coagulation. The tailored treatment approach and 'smooth motor' tip delivery, reduce discomfort and minimize damage to the skin's surface for optimal patient experience and outcomes.

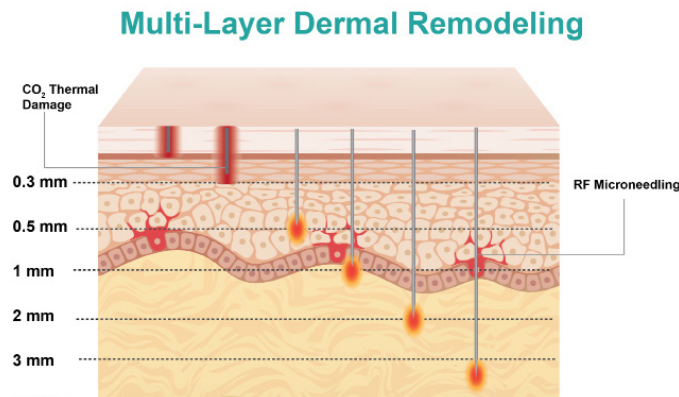

Deep Dermal Remodeling on All Skin Types


With each pulse, non-insulated or semi-insulated needles with 64 or 25 gold-plated pins, deliver RF energy to varying depths in the skin, adjustable in 0.1 mm increments from 0.5 mm to 3.5 mm.

CO₂ Fractional Laser Resurfacing

Secret PRO features a high-powered fractional CO₂ laser scanner that delivers consistent, uniform energy in short pulses. The ablative thermal effect resurfaces the skin and regenerates collagen for dermal remodeling.

Patterns and Modes


With one of the most dynamic scanners on the market, Secret PRO has a 100 micrometer beam size with 6 scan patterns and 9 scan modes to more easily treat difficult to target areas.

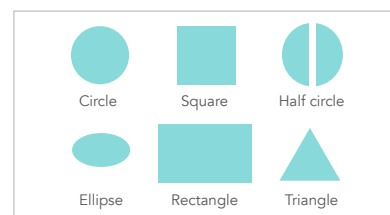

Treatment Area: 0.01 - 4 cm²
Dot Distance: 0.1 - 2.0 mm

i-Stack

Stack 1-10 laser pulses for adjustable depth with the i-Stack function to spot treat difficult scars¹ or deliver a customized full-face treatment, tailored to your patient's needs and tolerance for downtime.

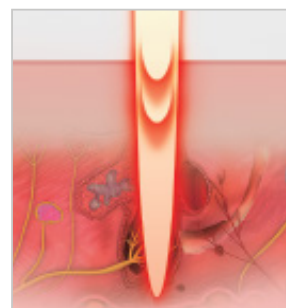

Four Unique Treatments

Address your patient's unique needs with the four treatment options below.



Refine

Radio Frequency Microneedling

A series of customizable RF microneedling treatments to provide deep dermal remodeling. Energy is precisely delivered at varying depths to create hemostasis and coagulation. This promotes collagen remodeling and elastin growth to address common skin concerns such as the signs of aging, scars, loss of radiance, skin tone, quality, and stretch marks.^{2, 3, 4}



Ultra-Light

CO₂ Light Laser Peel

A fractionally ablative, micro-peel for patients seeking a renewed, refreshed look, but don't need extensive CO₂ laser energy. The treatment with a shorter downtime targets photoaging, uneven texture, scars, and skin quality to deliver results in a single treatment.^{1, 5, 6}



Ultra

CO₂ Fractional Resurfacing Laser

For patients who desire a more aggressive treatment, this procedure uses CO₂ laser energy to target photoaging, skin quality, scars, and deliver results in a single treatment.^{1, 5, 6}



Revive

Radio Frequency Microneedling & CO₂ Laser

Reverse the signs of aging with the Revive treatment that layers RF microneedling and CO₂ laser energy in a single session. This procedure stimulates collagen, hyaluronic acid production, and improves overall skin quality.^{3, 6}

Secret[™] PRO

Real Results



Before
Somenek + Pittman MD



After
1 Treatment



Before
Somenek + Pittman MD



After
1 Treatment



Before
A. Bhatia, M.D.



After
1 Treatment, Subcision, and Fillers



Before
T. Lapa, M.D. & R. Mandavia, M.D.



After
1 Treatment



Before
B. Patel, M.D.



After
1 Treatment



Before
M. Aggarwal, M.D.



After
1 Treatment



Before
M. Aggarwal, M.D.



After
1 Treatment



Before
Somenek + Pittman MD



After
1 Treatment

Secret™ PRO

Technical Specifications

Technology	CO ₂ Fractional Laser	RF Microneedling
Source	10,600 nm wavelength laser	2 MHz bipolar radio frequency
Power	Up to 30 W (up to 175 mJ)	25 W
Pulse Duration	20 - 5,000 microseconds	50 - 950 milliseconds at 0 - 100% intensity
Repetition	0.2 / 0.5 / 0.8 / 1 sec / 2 sec / single	0.2 / 0.5 / 0.8 / 1 sec / 2 sec / single
Treatment Area	0.01 - 4 cm ²	25 pin 1 cm ² / 64 pin 2.25 cm ²
Spot Size	Fractional: 100 - 120 microns Surgical: Min. 0.2 - max. 1.3 mm	
Distance	0.1 - 2.0 mm	
Handpiece Options	Fractional and surgical available	1 handpiece fits all tip sizes
Adjustable Depths	Stack 1 - 10 laser pulses for adjustable depth with i-Stack function	0.5 - 3.5 mm (adjustable in 0.1 mm increments)
Needle Tips		25 non-insulated / 25 semi-insulated / 64 semi-insulated
Needle Thinness		0.25 mm gold plated

About Cutera

Cutera is a leading provider of laser and other energy-based aesthetic systems for practitioners worldwide. Since 1998, Cutera has been developing innovative, easy-to-use products that enable physicians and other qualified practitioners to offer safe and effective aesthetic treatments to their patients.



Cutera, Inc. World Headquarters
3240 Bayshore Boulevard, Brisbane, California 94005 USA
888-4-CUTERA or 415-657-5500 | info@cutera.com | cutera.com | Nasdaq CUTR

- ¹ Magnani, Lauren Rose, and Eric S. Schweiger. "Fractional CO₂ Lasers for the Treatment of Atrophic Acne Scars: A Review of the Literature." *Journal of Cosmetic and Laser Therapy*, vol. 16, no. 2, 2013, pp. 48-56., <https://doi.org/10.3109/14764172.2013.854639>.
- ² Juhasz, Margit L.W., and Joel L. Cohen. "Microneedling for the Treatment of Scars: An Update for Clinicians." *Clinical, Cosmetic and Investigational Dermatology*, Volume 13, 2020, pp. 997-1003., <https://doi.org/10.2147/ccid.s267192>.
- ³ Yadav, Savita, and Aashim Singh. "Microneedling: Advances and Widening Horizons." *Indian Dermatology Online Journal*, vol. 7, no. 4, 2016, p. 244., <https://doi.org/10.4103/2229-5178.185468>.
- ⁴ Cho, Soo Ick, et al. "Evaluation of the Clinical Efficacy of Fractional Radiofrequency Microneedle Treatment in Acne Scars and Large Facial Pores." *Dermatologic Surgery*, vol. 38, no. 7, 2012, pp. 1017-1024., <https://doi.org/10.1111/j.1524-4725.2012.02402.x>.
- ⁵ Ryu, Han-Won, et al. "Clinical Improvement of Striae Distensae in Korean Patients Using a Combination of Fractionated Microneedle Radiofrequency and Fractional Carbon Dioxide Laser." *Dermatologic Surgery*, vol. 39, no. 10, 2013, pp. 1452-1458., <https://doi.org/10.1111/dsu.12268>.
- ⁶ Tierney, Emily, and William Hanke. "Fractionated Carbon Dioxide Laser Treatment of Photoaging: Prospective Study in 45 Patients and Review of the Literature." *Dermatologic Surgery*, vol. 37, no. 9, 2011, pp. 1279-1290., <https://doi.org/10.1111/j.1524-4725.2011.02082.x>.

©2022 Cutera, Inc. All rights reserved. Model; Not actual patient. AP003543 rD (02/22)

